

medel[®]

FAMILY

AEROSOL THERAPY SYSTEM



Efficient
Low noise



Compact design



Easy to use




5-year warranty

Medel **FAMILY** is a piston-type compressor nebulizer system.


The convenient and functional design makes Medel FAMILY the ideal device for the whole family. The nebulizer ensures maximum speed and effectiveness; it comes with a convenient bag which allows storing all the accessories after use. It is equipped with a complete range of accessories.

Medical Device Class IIa.

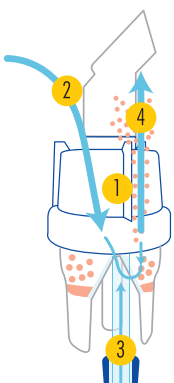

Medeljet Evo

Thanks to the Venturi effect, the air **2** is taken from the environment through top opening **1** in addition to the air supplied by the compressor **3**, thereby increasing nebulization speed and decreasing treatment time. **4**

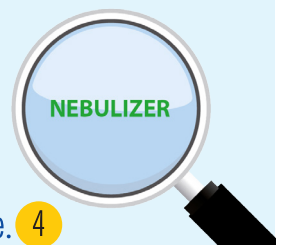

FAMILY

Medel **FAMILY** is supplied with a complete nebulization kit:
Medeljet Evo nebulizer, mouthpiece, nosepiece,
adult & pediatric masks, spare filters and air tube.



Technical specifications	FAMILY
Flow	5,3 l/min
Aerosol output	0,171 ml
Aerosol output rate	0,082 ml/min
Nebulization rate	0,4 ml/min
Particle size (MMAD)	3,16 µm
Fine Particle Fraction (FPF)	71,20%
Compressor maximum pressure	2,2 bar
Noise level	52dBA
Minimum filling volume	2 ml
Maximum filling volume	8 ml
Operating mode	30 min. on / 30 min. off
Weight	1,5 kg
Size	188 mm × Ø 136 mm

CODE 92459 Medel FAMILY

CE 0123



www.medel.it